

LUTRONIC

INFINI

High Intensity Focused RF

REVERSE TIME WITH
3D VOLUMIZATION



Now You Can Safely Treat All Skin Types... Year Round

HIGH INTENSITY FOCUSED RF

The LUTRONIC INFINI is the high intensity focused RF solution to create targeted and robust coagulation zones in aging, sagging skin with pigmentation, wrinkles and laxity. Complete control over treatment depth and energy levels enables a multi-layered targeted approach for truly customized treatments, regardless of patient's skin type or tanned skin. Now you can provide the best solution for results your patients demand year round.

UNIQUE TECHNOLOGY, AMAZING RESULTS

Infini's proven and impressive design delivers enhanced wrinkle reduction through 3-dimensional energy delivery to volumize tissue. High intensity focused RF delivered via insulated gold-coated microneedles creates precise and controllable fractionated coagulation zones within a specific layer of dermis. Adjustable depth control allows for customized and reproducible treatments of delicate areas. Infini effectively lowers the risk of PIH in dark skin types and can be safely used on patients with melasma and other types of pigmented lesions. Lutronic's patented ComfortFlow™ technology maximizes comfort and patient satisfaction.


RESULTS YOU CAN COUNT ON



Courtesy of Matteo Tretti Clementoni, MD, Italy



Courtesy of Steve Weiner, MD, USA

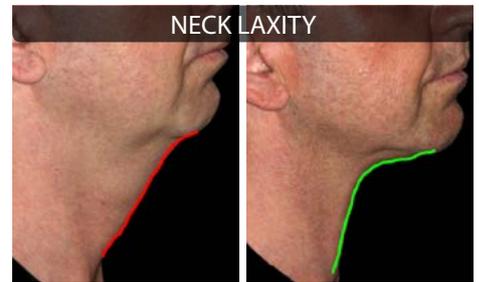

Courtesy of Matteo Tretti Clementoni, MD, Italy



Courtesy of Ivan Perez Haded, MD, Colombia



Courtesy of S. S. Savant, MD, India



Courtesy of J. Sakhiya, MD, India

MAJOR INDICATIONS

- Skin Tightening
- Periorbital Wrinkles
- Skin Rejuvenation
- Laxity in Jowl and Neck
- Active Inflammatory Acne
- Enlarged Pores
- Acne Scar
- Excessive Sebum Secretion
- Hyperhidrosis/Osmidrosis
- Stretch Mark

WHY INFINI?

- Proven; more than 600 devices sold globally
- Creates thousands of targeted coagulation zones per treatment to smooth wrinkles and improve aging, sagging skin
- Less pain and cost than micro-focused ultrasound
- More effective with broader range of applications than other non-ablative laser treatments
- Low cost consumables for maximized ROI
- Fast treatment; deliver more than 50k Thermal Coagulation Points (TCPs) in one 45 minute session
- Quick and easy recovery, less downtime than CO2
- 2 in 1 expandable platform: MFR & SFR
- Featured in ELLE magazine as the 9 Game-Changing Beauty Treatments that Defined 2015


WHAT PHYSICIANS ARE SAYING

"Finally! A Radiofrequency System That Makes Sense: The Infini... Infini has the right combination of the transmission of energy through the insulated microneedle tips and variable depths taken together with the flexible balance between power levels and exposure times, making Infini the ideal 'RF Skin Tightener'". Steve Weiner, MD, USA


INNOVATIVE COMFORTFLOW TECHNOLOGY



Insulated gold-coated microneedles

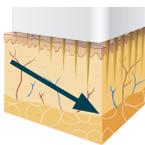

Built-in ComfortFlow technology enables smoother penetration via sequential microneedle insertion which minimizes patient discomfort for quick and easy recovery.


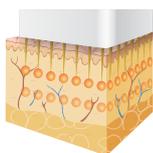


Adjustable depth control allows for multiple passes and treatment on delicate areas.

UNIQUE DUAL CHANNELS FOR SFR™






RF delivery with gentle skin contact

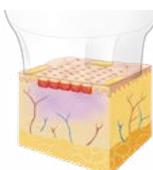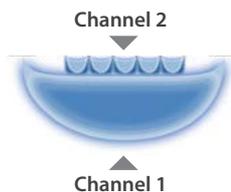

Unique dual channels of SFR noninvasively delivers RF energy to the epidermis and superficial dermis with synergistic effect, resulting in luminous skin and controlled tightening.

IMPROVED OUTCOMES

Infini's consistent, even thermal distribution ensures more consistent results and improved outcomes.

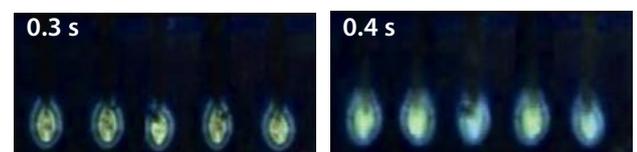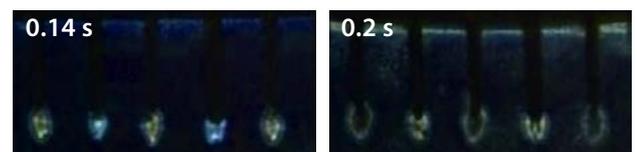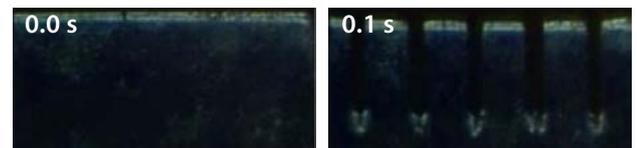


Microneedle tracks and thermal distribution of RF energy in the skin phantom model for the 3.5 mm depth at power level 8 and with different exposure times as shown (in seconds). Images were captured immediately after withdrawal of the microneedles.

Courtesy of Rudolf Verdaasdonk, PhD, et al., Fysica en Medische Technologie (FMT) Vrije University, Amsterdam, Netherlands

PROVEN CLINICALLY EFFECTIVE

Infini is clinically proven, with studies to back the technology. Regardless of depth you can be confident that you are creating consistent coagulation zones which provide outcomes you can count on.

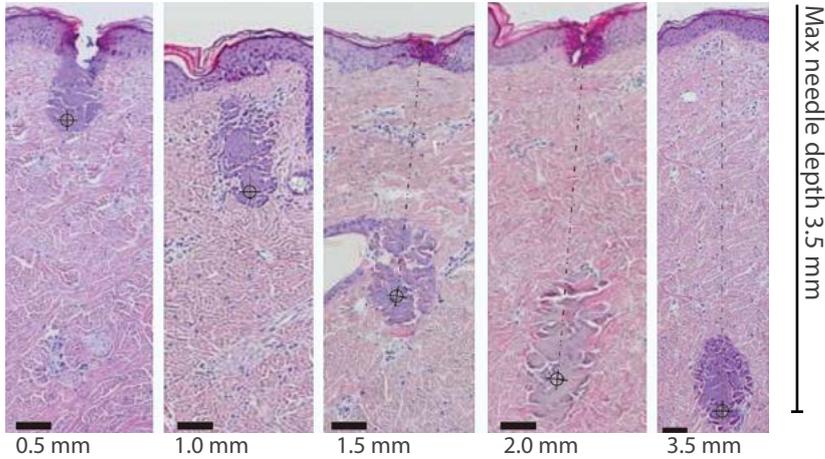

0.5 mm 1.0 mm 1.5 mm 2.0 mm 3.5 mm

All depths treated at Level 4 for 100 ms

Courtesy of SB Cho, MD, Dermatologist, Yonsei University, S. Korea

BEHIND HIGH INTENSITY FOCUSED RF

The robust Infini has a long list of benefits that support you and your patients' needs.

- Ergonomically designed lightweight handpiece
- Precise targeting with adjustable microneedles over 3.5 mm
- Foot or finger switch operated
- Intuitive touch screen GUI
- Patented ComfortFlow™ technology
- Expandable platform with two revenue streams

SYSTEM SPECIFICATIONS

System Type	Bipolar RF (Radiofrequency)
Frequency	1 MHz
Total RF Treatment Level	Up to Level 20
Operation Mode	Microneedle Fractional RF (MFR) Superficial Fractional RF (SFR)
Treatment Duration (Time)	10 ms - 1,000 ms
Tips	MFR : Microneedle electrodes, up to 3.5 mm (Adjustable depth) SFR : Up to 144 pins
Display	Color Touch LCD (Total Shot Count Displayed)
Dual Switch	Foot & Finger Switch
Weight	28 kg
Dimension	362 mm (W) x 409 mm (L) x 1713 mm (H)



Lutronic Corporation

Lutronic Center, 219 Sowon-ro, Deogyang-gu, Goyang-si, Gyeonggi-do, KOREA
TEL. +82 31 908 3440 - FAX. +82 1600 3660 - office@lutronic.com

riyadh +966 11 465 0371
qassim +966 16 326 5404
kuwait +965 2 225 0206
beirut +961 1 998 189

jeddah +966 12 639 8200
khamis mushait +966 17 237 0072
manama +973 1733 1046
amman +962 6 553 9977

madina +966 14 813 2302
dubai +971 4 451 4455
doha +974 4441 6893
others countries +971 4 451 4455

khobar +966 13 894 5051
abu dhabi +971 2 674 6493
muscat +968 24 297 816