

New!
with **1064 nm**

Technology You Trust Innovation You Need



SPTL-1

First FDA-cleared PDL



Vbeam Classic

Proprietary
595 nm wavelength



Vbeam Perfecta

8 micropulses



Vbeam Prima

- Dual wavelength
- Dual cooling
- More user-friendly features
- Better treatment experience




Vbeam® Prima

A recognized leader in vascular treatments, Candela now brings you the latest pulsed-dye laser with more enhancements than ever before.

Clinically proven results across indications

Versatile treatments

Vbeam Prima is a highly efficient pulsed-dye laser (PDL) used by healthcare providers all over the world to treat a variety of indications for both face and body, including benign vascular, pigmented, and certain non-pigmented lesions, with high patient tolerability and a low incidence of side effects.¹⁻¹²

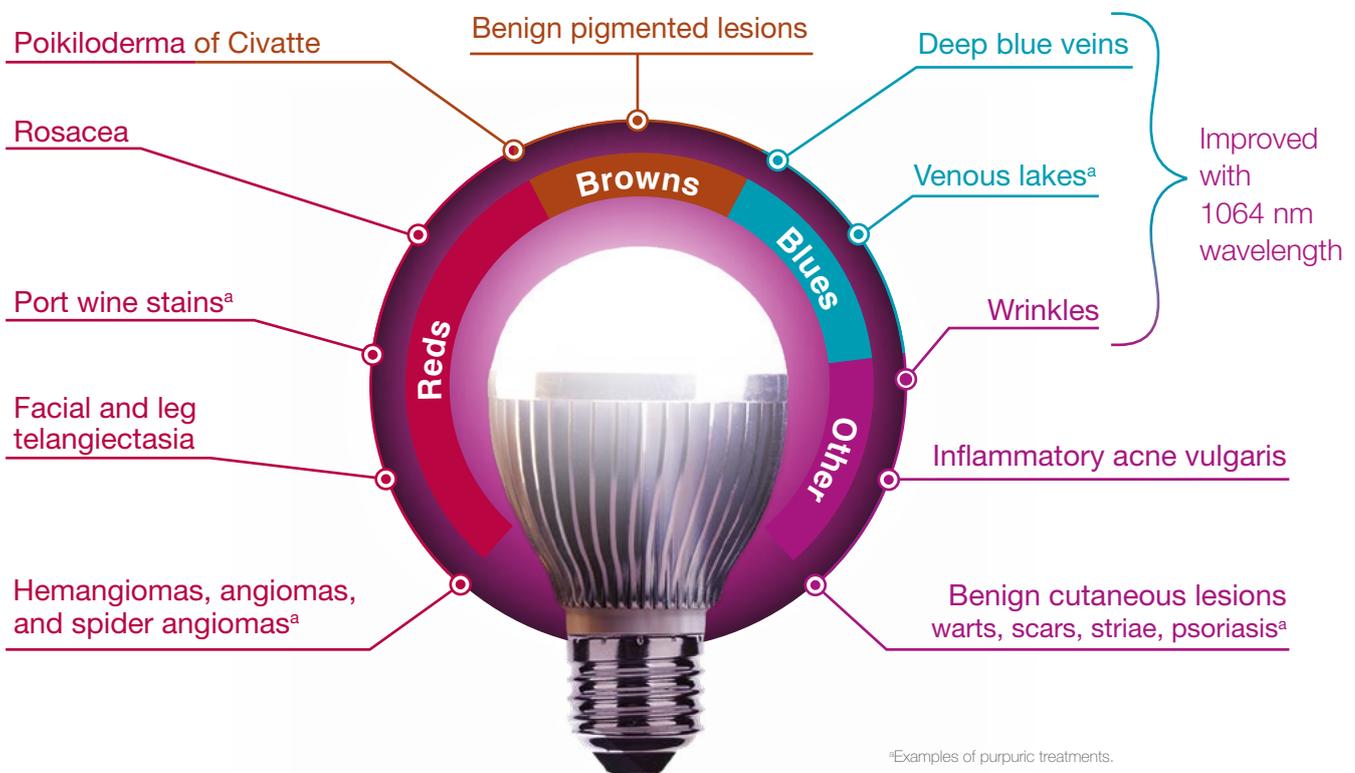


A legacy of innovation

<p>1983 Anderson/Parrish MGH/Harvard</p>	<p>1988 SPTL-1 577 nm, first FDA-cleared PDL</p>	<p>1994 SPTL-1b 585 nm</p>	<p>1996 ScleroPlus® 585 nm, 590 nm, 595 nm, 600 nm</p>
	<p><small>SPTL, selective photothermolysis.</small></p>		

Two wavelengths, more treatment parameters

Proprietary 595 nm wavelength


The 595 nm wavelength deeply penetrates the skin to reach targeted blood vessels.⁴ Its energy is absorbed by oxyhemoglobin to coagulate and clear vessels with more tolerability and fewer instances of melanin absorption.¹³

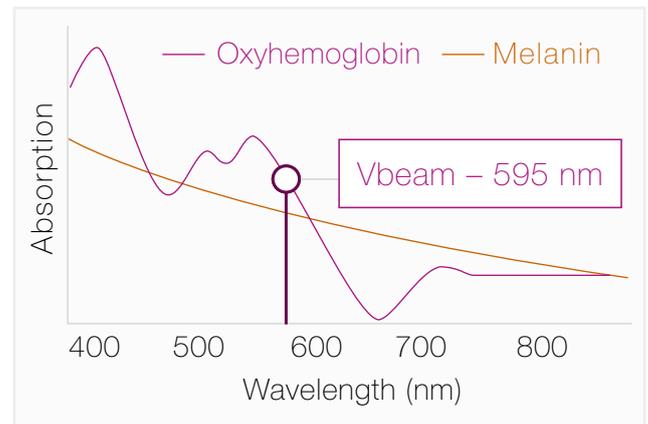

Additional 1064 nm wavelength

Treat deep blue veins, venous lakes, and wrinkles!¹²

- Treats blue veins across face and body
- Minimizes the appearance of wrinkles

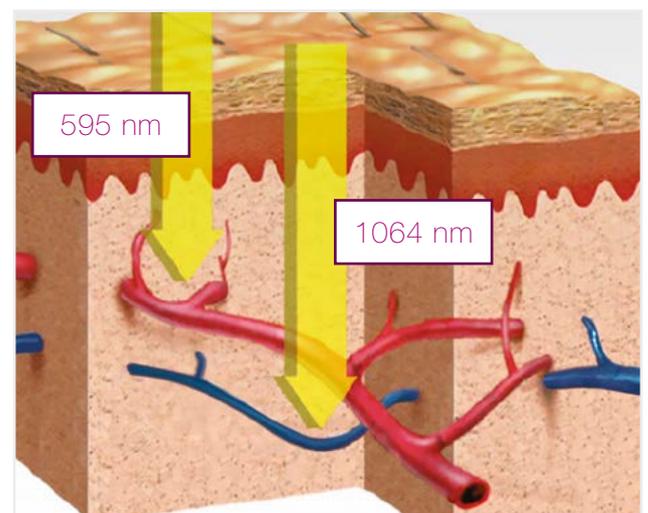

2000

Vbeam Classic
595 nm

2006

Vbeam Perfecta®
595 nm
8 micropulses


2018

Vbeam Prima
595 nm and 1064 nm



New features that benefit both the provider and patient experience

Achieve greater results in less time

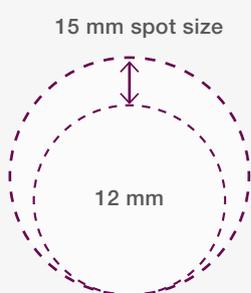


Larger max treatment spot size and greater max energy with 595 nm

- Faster, more efficient treatments
- A treatment that requires 100 pulses with a 12 mm spot size can require as few as 64 pulses with a 15 mm spot size

	VBEAM PERFECTA	VBEAM PRIMA
Maximum energy	8 J	12 J
Maximum spot size	12 mm	15 mm
Area covered	1.13 cm ²	1.77 cm²
Fluence at 12 mm	7.1 J/cm ²	10.6 J/cm²




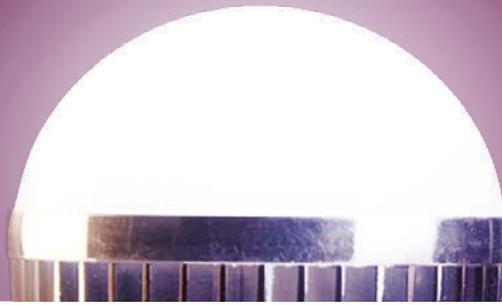

EverCOOL

Patented cryogen-based Dynamic Cooling Device™ (DCD™) and EverCool™ contact cooling options

- Two types of cooling for maximum versatility and epidermal protection
- Patented DCD scales with fluence to automatically administer consistent epidermal protection
- With EverCool, you can treat vascular and pigmented lesions at the same time

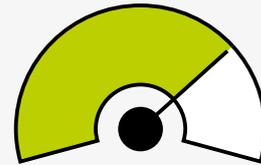

Smart dye life management

- Dye life meter provides data on remaining dye levels
- Avoid unexpected treatment interruptions from dye loss
- Prompts user to schedule preventative maintenance



Calibration on first use

Device is patient- and treatment-ready for the entire day, reducing downtime between patients; no need to recalibrate when changing spot sizes or fluences



Wi-Fi connectivity

More accurate remote service diagnostics for faster equipment servicing



Zoom handpiece

A versatile handheld applicator that allows for very targeted spot size adjustments in increments as small as 0.5 mm




Guided user interface


Intuitive, easy-to-use software with quick access to saved, favorite treatment settings

Powerful results begin with Vbeam Prima


Vbeam Prima is the trusted, proven PDL across indications¹⁴

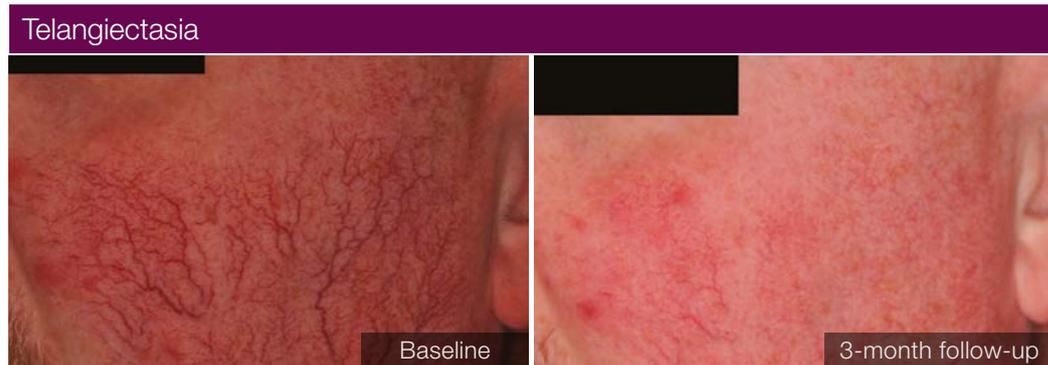

Photos are unretouched. Patient treated with Vbeam Prima; individual results may vary. Photos courtesy of E. Victor Ross, MD.

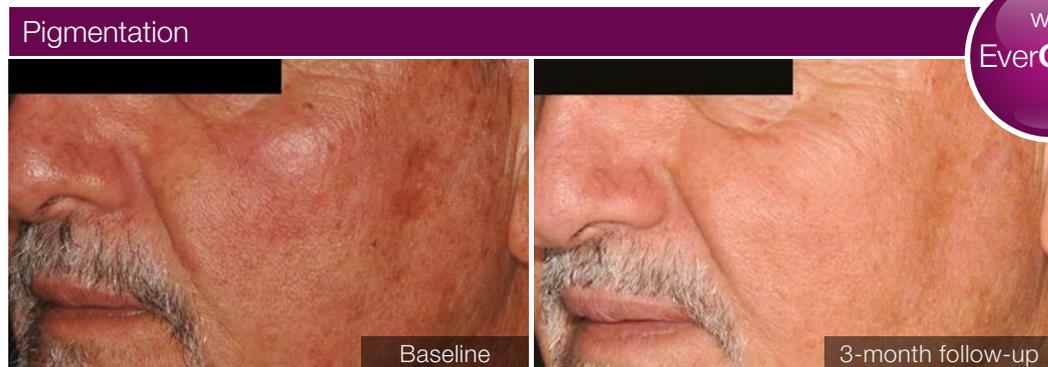


Photos are unretouched. Patient treated with Vbeam Prima; individual results may vary. Photos courtesy of E. Victor Ross, MD.

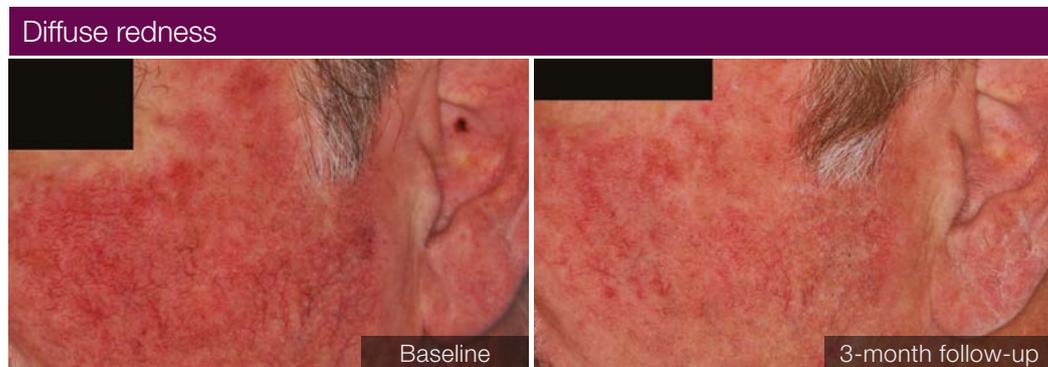

Photos are unretouched. Patient treated with Vbeam Prima; individual results may vary. Photos courtesy of E. Victor Ross, MD.

Demand the best for your patients

Rosacea affects an estimated 415 million people worldwide¹⁵

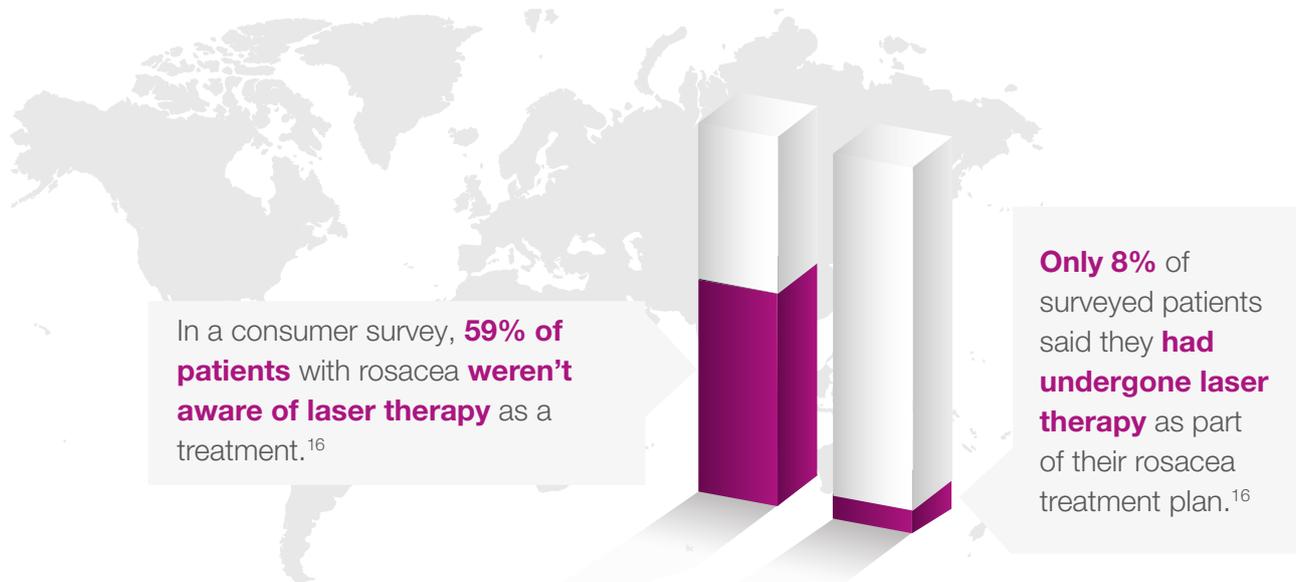

Significant improvements in rosacea symptoms with Vbeam Prima¹⁷

In a clinical study, rosacea symptoms were significantly reduced from moderate to mild after 4 treatments.

- ☑ Difficult-to-treat erythema was reduced from severe to mild
- ☑ ~90% of patients had >40% rosacea improvement



of patients reported being very satisfied with results¹⁷



Photos are unretouched. Patient treated with Vbeam Prima; individual results may vary. Photos courtesy of Eric F. Bernstein, MD, MSE.

Minimal adverse side effects reported¹⁷

On a scale of 1 to 10 (10=max), patients reported an average pain score of 5.6 ± 1.8 . Most patients reported side effects such as mild edema, mild to moderate erythema, and mild to moderate bruising. These effects resolved within a few days after treatment without intervention.

Experience the benefits of Vbeam Prima for versatile treatments with proven results

- **Larger 15 mm spot size and 50% greater power at 12 J with 595 nm**— for faster, more efficient treatments
- **Calibration on first use** — for less downtime between treatments, and more efficient treatments when using a variety of settings
- **Smart dye life management**— for real-time data on remaining dye levels
- **1064 nm wavelength**— for deep blue vessel clearance and wrinkle reduction
- **Patented cryogen-based DCD and EverCool contact cooling**— for maximum epidermal protection and treatment versatility

PRODUCT SPECIFICATIONS

Laser wavelengths	<ul style="list-style-type: none"> › PDL: 595 nm › Nd:YAG laser: 1064 nm
Laser pulse repetition rate	<ul style="list-style-type: none"> › PDL: Up to 1.5 Hz › Nd:YAG laser: Up to 10 Hz
Laser pulse duration	<ul style="list-style-type: none"> › PDL: 0.45-40 ms › Nd:YAG laser: 0.5-40 ms
Maximum pulse energy	<ul style="list-style-type: none"> › PDL: 12 J › Nd:YAG laser: 30 J
Method of optical output	Lens-coupled optical fiber with user-selectable spot sizes
Networking method	Wi-Fi
Dimensions (H x W x D)	53 x 20 x 33 in / 135 x 51 x 84 cm
Weight	280 lbs / 127 kg
Electrical requirements	220-240 VAC, 24 A max, 50/60 Hz, single phase
Cooling methods	Patented Dynamic Cooling Device (DCD) and EverCool contact cooling

Patented Dynamic Cooling Device (DCD)

Integrated controls, cryogen container, and handpiece with distance gauge

Cryogen	HFC 134a
DCD spray duration	User-adjustable range: 20-100 ms
DCD delay duration	User-adjustable range: 10-100 ms
DCD post-spray duration	User-adjustable range: 0-50 ms
Beam spot sizes	<ul style="list-style-type: none"> › 3-15 mm Zoom handpiece › 3 x 10 mm

Patented EverCool contact cooling

User-controlled, adjustable sapphire cooling tip for use pre, during, and post pulse

Temperature range	10°C-20°C
Beam spot size	3-15 mm Zoom handpiece

For more information about how Vbeam Prima* may help you achieve your practice goals, visit [syneron-candela.com](https://www.syneron-candela.com)

1. Vbeam 510(k) clearance (K033461), January 2004. 2. Vbeam 510(k) clearance for pigmented lesion handpiece accessory (K051359), July 2005. 3. Bernstein EF, Kligman A. Rosacea treatment using the new-generation, high-energy, 595 nm, long pulse-duration pulsed-dye laser. *Lasers Surg Med.* 2008;40(4):233-239. 4. Woo SH, Ahn HH, Kim SN, Kye YC. Treatment of vascular skin lesions with the variable-pulse 595 nm pulsed dye laser. *Dermatol Surg.* 2006;32(1):41-48. 5. Chapas AM, Eickhorst K, Geronemus RG. Efficacy of early treatment of facial port wine stains in newborns: a review of 49 cases. *Lasers Surg Med.* 2007;39(7):563-568. 6. Jasim ZF, Woo WK, Handley JM. Long-pulsed (6-ms) pulsed dye laser treatment of rosacea-associated telangiectasia using subpurpuric clinical threshold. *Dermatol Surg.* 2004;30(1):37-40. 7. Jørgensen GF, Hedelund L, Hæderdal M. Long-pulsed dye laser versus intense pulsed light for photodamaged skin: a randomized split-face trial with blinded response evaluation. *Lasers Surg Med.* 2008;40(5):293-299. 8. Halachmi S, Israeli H, Ben-Amital D, Lapidot M. Treatment of the skin manifestations of hereditary hemorrhagic telangiectasia with pulsed dye laser. *Lasers Med Sci.* 2014;29(1):321-324. 9. Yu W, Ma G, Qiu Y, et al. Prospective comparison treatment of 595-nm pulsed-dye lasers for virgin port-wine stain. *Br J Dermatol.* 2015;172(3):684-691. 10. Galeckas KJ, Ross EV, Uebelhoer NS. A pulsed dye laser with a 10-mm beam diameter and a pigmented lesion window for purpura-free photorejuvenation. *Dermatol Surg.* 2007;34:1-6. 11. Madan V, Ferguson J. Using the ultra-long pulse width pulsed dye laser and elliptical spot to treat resistant nasal telangiectasia. *Lasers Med Sci.* 2010;25(1):151-154. 12. Indications for 1064 nm wavelength. Candela, data on file. 13. Bernstein EF. The pulsed-dye laser for treatment of cutaneous conditions. *G Ital Dermatol Venerol.* 2009;144(5):557-572. 14. Ross EV. Vbeam Prima before and after photos. Candela, data on file. 15. National Rosacea Society website. <https://www.rosacea.org/weblog/415-million-people-affected-rosacea-worldwide>. Accessed May 8, 2018. 16. Consumer rosacea laser attitude and behavior exploratory: final report. September 7, 2017. BuzzBack Market Research. 17. Bernstein EF, Schomacker K, Paranjape A, Jones CJ. Pulsed dye laser treatment of rosacea using a novel 15 mm spot size. Candela, data on file.



© 2019 Candela Corporation. This material contains registered trademarks, trade-names and brand names of Candela Corporation and its affiliates, including Syneron, Candela, DCD, Dynamic Cooling Device, EverCool, ScleroPlus, Vbeam, and Vbeam Perfecta. All other trademarks are the property of their respective owners. All rights reserved. PB98813EN

riyyadh +966 11 465 0371 jeddah +966 12 639 8200 madina +966 14 813 2302 khobar +966 13 894 5051
 qassim +966 16 326 5404 khamis mushait +966 17 237 0072 dubai +971 4 451 4455 abu dhabi +971 2 674 6493
 kuwait +965 2 225 0206 manama +973 1733 1046 doha +974 4441 6893 muscat +968 24 297 816
 beirut +961 1 998 189 amman +962 6 553 9977 others countries +971 4 451 4455

imdad | امداد | Qualified Solution