



# Guide to Tips & Handpieces

## Tips: Overall Naming Convention

**Tip Name** (what each character means)

**Possible Values for Each Position**

1st character = size of treatment area	G = General use	M = Macro size	P = Precision
2nd & 3rd characters = number of pins	14 = rectangular 2 rows of 7 pins	16 = square 4 rows of 4 pins	49 = square 7 rows of 7 pins
4th character = target depth	D = Dermal		

## Common Applications and Handpiece Compatibility

Type of Tip	Coverage Area	Common Applications	Standard (BLUE) Handpiece	Standard (RED) Handpiece
G = General Use	10 x 10mm	Face, neck, body	YES	YES
M = Macro size	20 x 20mm	Off-face, larger areas	NO	YES
P = Precision	<i>Less than</i> 10 x 10mm	Periorbital, scars, striae, etc.	YES	YES



## Common Tips are Color-Coded

Tip Name	Package Color	Description of Tip (decoded name)	Standard HP	Macro HP
G49D	RED	General use (10x10mm), 49 pins (7x7 square), targeting dermis	YES	YES
M49D	GREEN	Macro size (20x20mm), 49 pins (7x7 square), targeting dermis	YES	YES
P14D	PURPLE	Precision size (less than 10x10mm), 14 pins (2 x 7 rectangle), targeting dermis	YES	YES
P16D	TEAL	Precision size (less than 10x10mm), 16 pins (4x4 square), targeting dermis	YES	YES

# Standard Handpiece



# Macro Handpiece

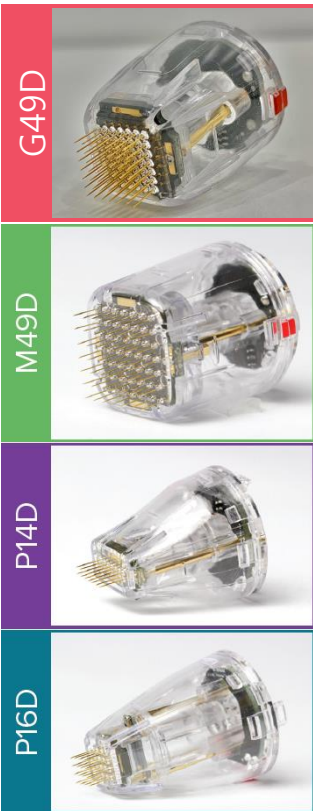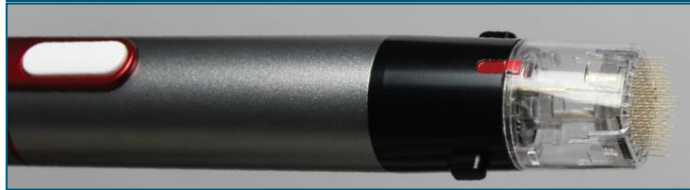

G49D	<b>LUTRONIC</b> Intelligent Care™
M49D	<b>LUTRONIC</b> Intelligent Care™
P14D	<b>LUTRONIC</b> Intelligent Care™
P16D	<b>LUTRONIC</b> Intelligent Care™